

lindkvist 



ULTRAVIEW WL



- ✓ Up to 9 operating hours without battery recharge
- ✓ Developed in cooperation with experienced surgeons
- ✓ Operating theatre approved light intensity.




The ULTRAVIEW WL headlight provides the surgeon with the ultimate in comfort, visibility and flexibility. It weighs in at only 310g and has a padded leather headband, which provides exceptional comfort. Each headlight comes with two spare batteries and a charger. The charge will last for up to nine hours, depending on the brightness setting.

ULTRAVIEW's key benefits to surgeons are:

- **Cordless technology**
- **Excellent visibility**
- **Lightweight and comfortable**
- **Three levels of brightness**
- **Consistent white light**
- **Very competitively priced**

ULTRAVIEW WL - Cable-free headlight innovation

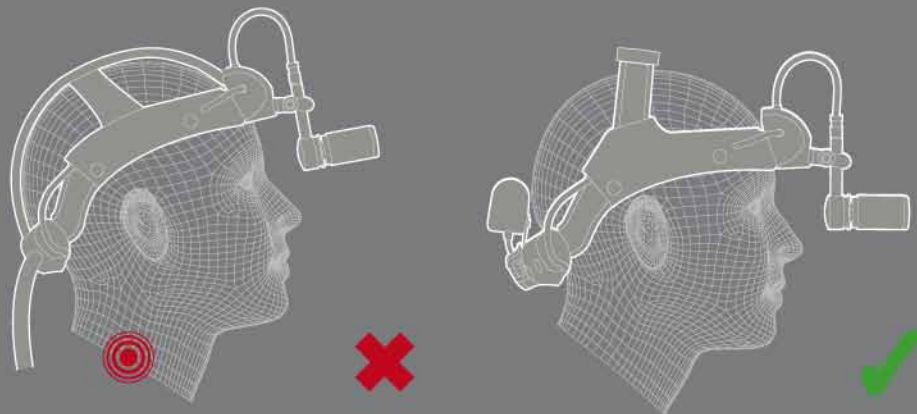

The ULTRAVIEW WL helps to eliminate neck pain caused by the weight of the leads in a conventional headlight and also eliminates the safety issues caused by trailing leads.



Key features and benefits



Easy to change Li-ion cell battery
Charging time 2.5 hours.



Control panel
Easy to use,
illuminated panel.



Push and twist adjustable headstrap
Fully adjustable by
pushing and turning.



Lightweight
The whole unit weighs only 310g.



Adjustable headband

To ensure comfort and stability.

**Water and liquid resistant
cable connections**

**Integrated ventilator
prevents overheating**

Padded leather headstrap

Ergonomically designed for comfort throughout extended periods of use.

CE marked to ISO 9001 13485



Fully adjustable light positioning

Easy to direct light beam.

Magnifying lenses

Clip-on magnifying lenses, +2D: 1.4x as standard other magnification strengths are available on request.



Focusable LED light

Twist barrel to zoom focus high power LED light with 3 levels of brightness.

Specification

Each **Standard-Set** includes 1 headlight, 2 batteries, 1 charger and 1 Aluminium Case

Ordering-no. 96100



ULTRAVIEW WL Headlight

Illumination	45,000Lux (at 25cm working distance) level 3 35,000Lux (at 25cm working distance) level 2 17,000Lux (at 25cm working distance) level 1
Colour temperature	White 6000 Kelvin
Illumination field	50mm to 75mm (at 25cm working distance)
Spot adjustment	Zooming control
Brightness	3 level control
LED life Span	Up to 50,000 hours
Operating time	9 hours (brightness level 1) 3.5 hours (brightness level 2) 2.5 hours (brightness level 3)
Output	3.7V
Weight	255g



ULTRAVIEW WL Rechargeable Batteries (2 included in standard set)

Ordering-no. 96120

Type	Li-ion (built-in protected circuit module)
Capacity	2400mA
Weight	55g

ULTRAVIEW WL Charger (incl. power supply unit)

Ordering.-no. 96130

Input	AC 100-240V, 50/60Hz
Output	4.2 V/1300mA
Charging time	120min


**ULTRAVIEW WL Magnifying Lenses**

Lens +2.0 diopters Ordering.-no. 96142

Lens +3.0 diopters Ordering.-no. 96143

Lens +4.0 diopters Ordering.-no. 96144



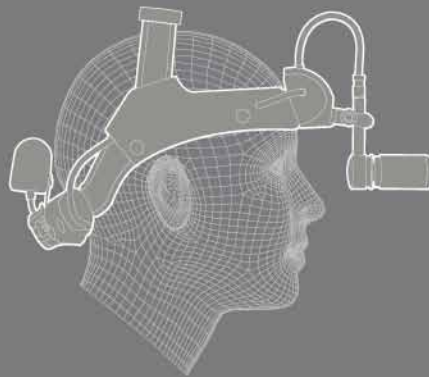

lindkvist 

nordiska 

Victoria Höhe 1 · D-51645 Gummersbach
Telefon +49(0)2261/9891-0 · Freecall 0800/8835000 · Telefax +49(0)2261/9891-111
info@nordiska.de · www.nordiska.de